

Ultra Shoulder Positioner[®]

For arthroscopic and open shoulder surgery



Enhanced Ease of Use

The **Ultra Shoulder Positioner**[®] (USP) provides a safe, convenient and stable method of maintaining patients in the “beach chair” position for arthroscopic open shoulder surgery and joint replacement. It is designed to easily attach to any general surgery table and suit the high demands of the operating room.

The USP provide many features that enhance its general ease of use. Its centrally located backrest ensures that unobstructed access is provided to the surgical site. The radiolucent backrest has a low contour for optimal shoulder access both anteriorly and posteriorly. All release levers and adjustment knobs are readily located for convenient access and easier positioning of the patient.



400 lbs.
weight capacity (181 kg)



The headrest wings on the USP retract for an obstruction-free space around the patient's head. The headrest is fully adjustable, easily accommodating anterior translation and lateral tilt.



⊕ Ultra Shoulder Positioner® Features

A. Adjustable Cervical Flexion

- Headrest accommodates kyphotic patient
- Headrest provides 4 inches (10.2 cm) of anterior translation

B. Headrest

- Headrest wings retract providing obstruction-free space for anesthetist
- Head strap and headrest wings provide secure head support
- Headrest tilts 30 degrees in both directions

C. Patient Care Kits


- Multiple, single-use patient care kits provide a variety of options including Shoulder Immobilizer and Rehab Shoulder Sling
- Multiflex™ Back Pad perforation provide tears to increase shoulder access on the operative side

D. Translating Arm Board and Lateral Support Pad Assembly

- Positioned arm support platform with translating arm board
- Lateral support pad provides lateral support and prevents patient shifting

E. Storage Cart

- Provides storage for the USP, accessories, and patient care kits
- Facilitates easy transport to and from surgery site

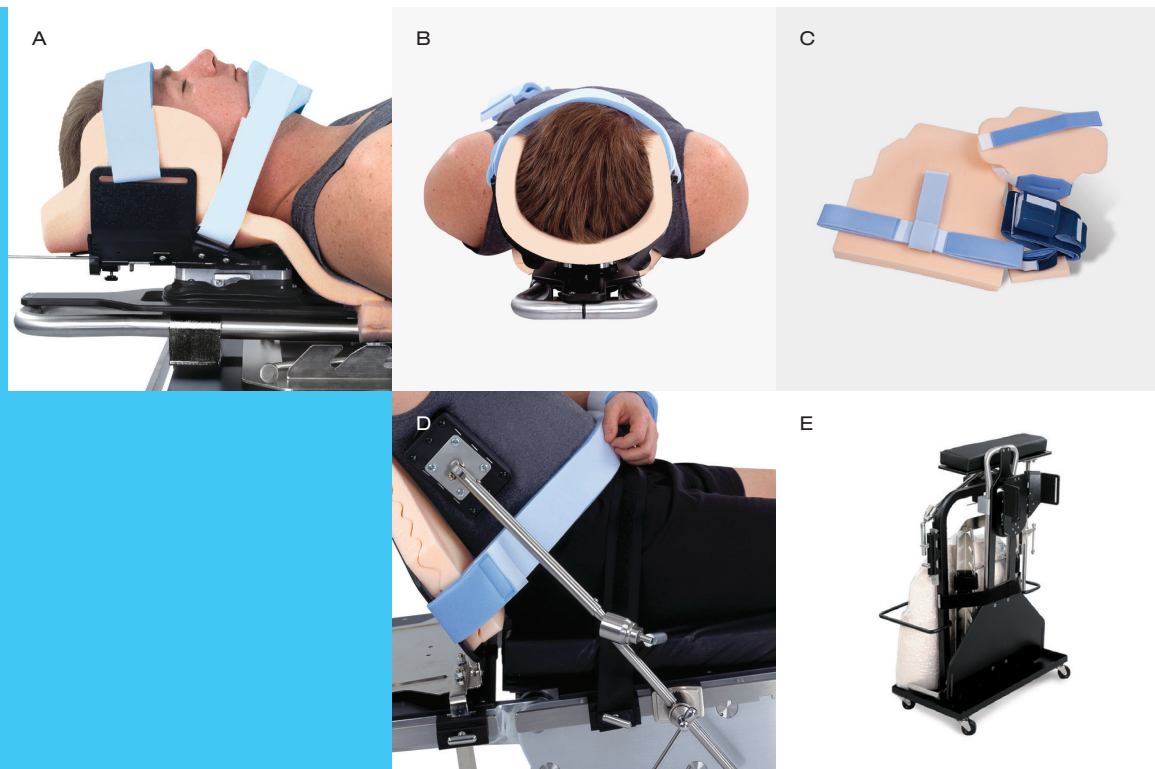

Specifications & Components



Specifications

Headrest Translation

Anterior Range	4.0 in. (10.2 cm)
Longitudinal Range	8.5 in. (21.5 cm)
Lateral Tilt	30 Degrees

Ultra Shoulder Positioner Weight	31 lbs. (14 kg)
----------------------------------	-----------------

Patient Sitting Height Range	30.8-39.2 in. (78.2 - 99.6 cm)
------------------------------	--------------------------------

Patient Weight Capacity	400 lbs. (181 kg)
-------------------------	-------------------

Ordering Information

REF

5370	Ultra Shoulder Positioner® and Cart
5358	Ultra Shoulder Positioner® with two Patient Care Kits (includes Rehab Sling)
5338-24	Ultra Shoulder Positioner® Cart
5338-100	Lateral Support Pad Assembly (includes Clark Socket)
5338-70	Translating Arm Board
5393	Clark Socket

Attaches to any general surgical table.

Individual Patient Care Kits

REF


5345	Patient Care Kit (3/cs)
5346	Patient Care Kit with Shoulder Immobilizer (3/cs)
5348	Patient Care Kit with Rehab Sling (3/cs)



30031 Ahern Avenue
Union City, CA 94587-1234 USA
Telephone: 510-429-1500
Toll Free: 800-777-4674
Fax: 510-429-8500
Outside USA: +1-510-429-1500
mizuho.com
newhipnews.com

EC REP

Emergo Europe
Prinsessegracht 20
2514 AP The Hague
The Netherlands

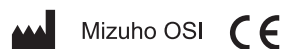

Note: Mizuho OSI is constantly improving its products. All specifications are subject to change without notice. Ultra Shoulder Positioner® is a registered trademark of Mizuho OSI. • Mizuho OSI is a Delaware Corporation. Manufactured in the USA.

2017 ©Mizuho OSI

REF NW0377 REV G